



TECHNICAL SPECIFICATIONS

Standard Seat Size: 17.5" wide x 15" deep (445 mm wide x 381 mm deep)

Larger Seat Option: 21.5" wide x 16" deep (546 mm wide x 406 mm deep)

Standard Footrest Size: 10" long x 9" wide (254 mm long x 229 mm wide)

Larger Footrest Option: 12" long x 11" wide (305 mm long x 279 mm wide)

Restore peace-of-mind and access in your home.

When you want the highest quality stairlift combined with a stylized appearance, you want the Bruno Elite. The Elite offers worry-free access to all levels of your home. And the Elite's stylized design blends with your home's environment. Consider custom upholstery the chair in your own fabric.

Let Bruno's Elite safely connect you to all levels of your home again.

- Made-in-the-USA quality
- Easy, reliable operation
- Smooth, quiet ride

Learn more at www.bruno.com.



ELITE INDOOR STRAIGHT | SRE 2010




Close installation to wall.



ELITE OPTIONS

Standard tan vinyl with five optional upholstery selections. Or, supply your own fabric and Bruno will custom upholster the chair for you.

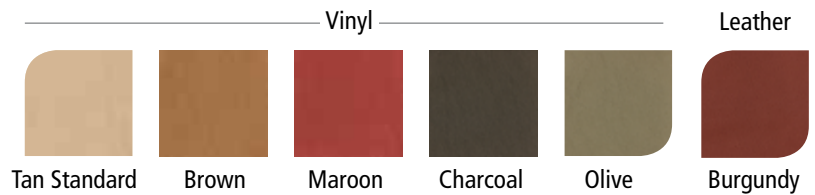

Learn more at www.bruno.com



Offset swivel seat makes it easy to get on/off.

ELITE FEATURES

- Weight capacity: 400 lb (181 kg)
- Plush, generous-size seat with multiple height adjustments
- Adjustable width armrests to 21.75" (552 mm) between arms
- Anodized rail with covered gear rack gives a sleek appearance
- Flip up arms, seat, and footrest create plenty of extra space on steps
- Adjustable footrest for maximum comfort
- Vertical rail design enables tight installation to within 6" (152mm) of wall
- Onboard audio/visual diagnostics for easy service
- Continuous charge strips along rail repower unit wherever it stops
- Two wireless call/sends allow unit to be controlled remotely
- Two 12v batteries ensure dependable performance even in a power outage
- Soft start/stop gives a smooth ride from start to finish
- Gold warranty: five-year major components; two-year parts

Luxury styling and comfort.



Power folding footrest automatically flips up/down when seat is raised/lowered.



Power or manual folding rails for narrow hallway or when doorway is at the bottom of stairs. Manual operation or push button automatic.



Power swivel seat for effortless exit. Controlled on chair arm or wireless call/send.



Larger seat and footrest for increased comfort.



About Bruno

Bruno Independent Living Aids has helped improve the lives of people with limited mobility for over 30 years. Engineered and manufactured in the USA, Bruno stairlifts and vertical platform lifts help people better access homes, public buildings and businesses, and Bruno's full line of scooter/powerchair lifts help people easily transport their mobility devices. In addition, Bruno's Valet Signature Seating is the world leader in revolutionary turning automotive seats to assist people in getting in and out of their vehicles more easily. A family-owned business, Bruno takes pride in their unwavering focus on quality and safety and their position as industry innovation leader.

Proudly made in the USA



Bruno Independent Living Aids, Inc.
1780 Executive Drive, Oconomowoc, WI 53066
www.bruno.com

© Bruno Independent Living Aids, Inc. 2015

MKT-SL-2